


pl 300

Cuscino in gel fluido Cushions in fluid gel

Codice Code	Misura (cm) Size (cm)
PL300/4040	40X40X4
PL300/4242	42X42X4
PL300/4545	45X45X4

*altre misure disponibili su richiesta
other sizes available on request*

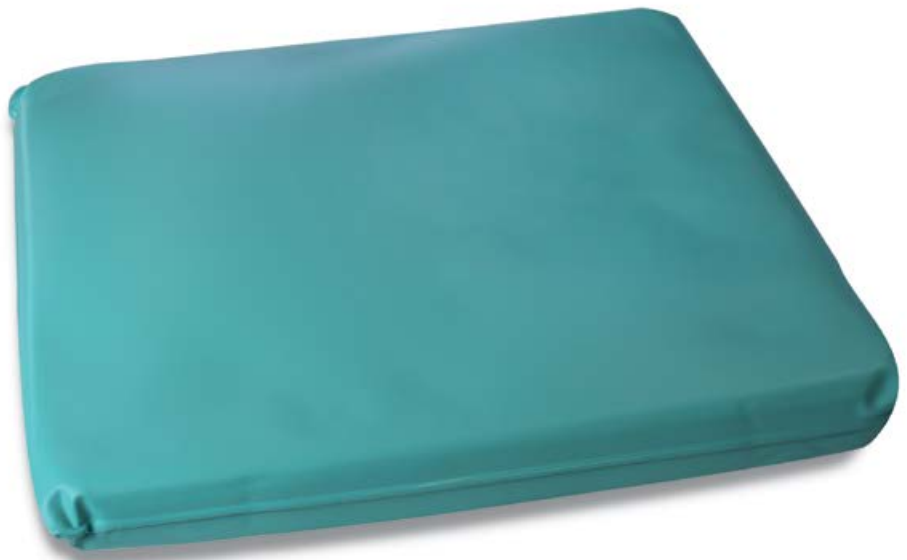

Cover

2 Bielastico
Bistretch

5 Bielastico + 3D + antiscivolo
Bistretch + 3D + antislip

Cuscino antidecubito costituito da gel fluido e da una base in poliuretano espanso, il tutto contenuto in una sacca in PVC.

La fluidità dello strato in gel superficiale consente un'ottima distribuzione del peso corporeo del paziente, diminuendo quindi la pressione sui capillari e prevenendone l'eventuale occlusione.

La base in poliuretano espanso permette al corpo di affondare nel cuscino, ma allo stesso tempo assicura un adeguato sostegno, così da evitare l'effetto bottom out e assicurare la stabilità della seduta.

Grado di rischio:

prevenzione.

Anti-decubitus cushion manufactured with fluid gel and a polyurethane foam base, those are contained in a hermetic PVC sack.

The gel surface layer fluidity permits an excellent patient's weight distribution, therefore decreasing the pressure on the capillaries, preventing their potential obstruction.

The polyurethane foam base permits to the body to sink in the cushion, but at the same time assures an appropriate support, thus to avoid the bottom-out effect and to assure the sitting stability.

Risk level:

prevention.

